


EUROTAPE

Connecting with care

Med/PV175N™

OR-tape for re-usable Drapes

100% Wash Soluble

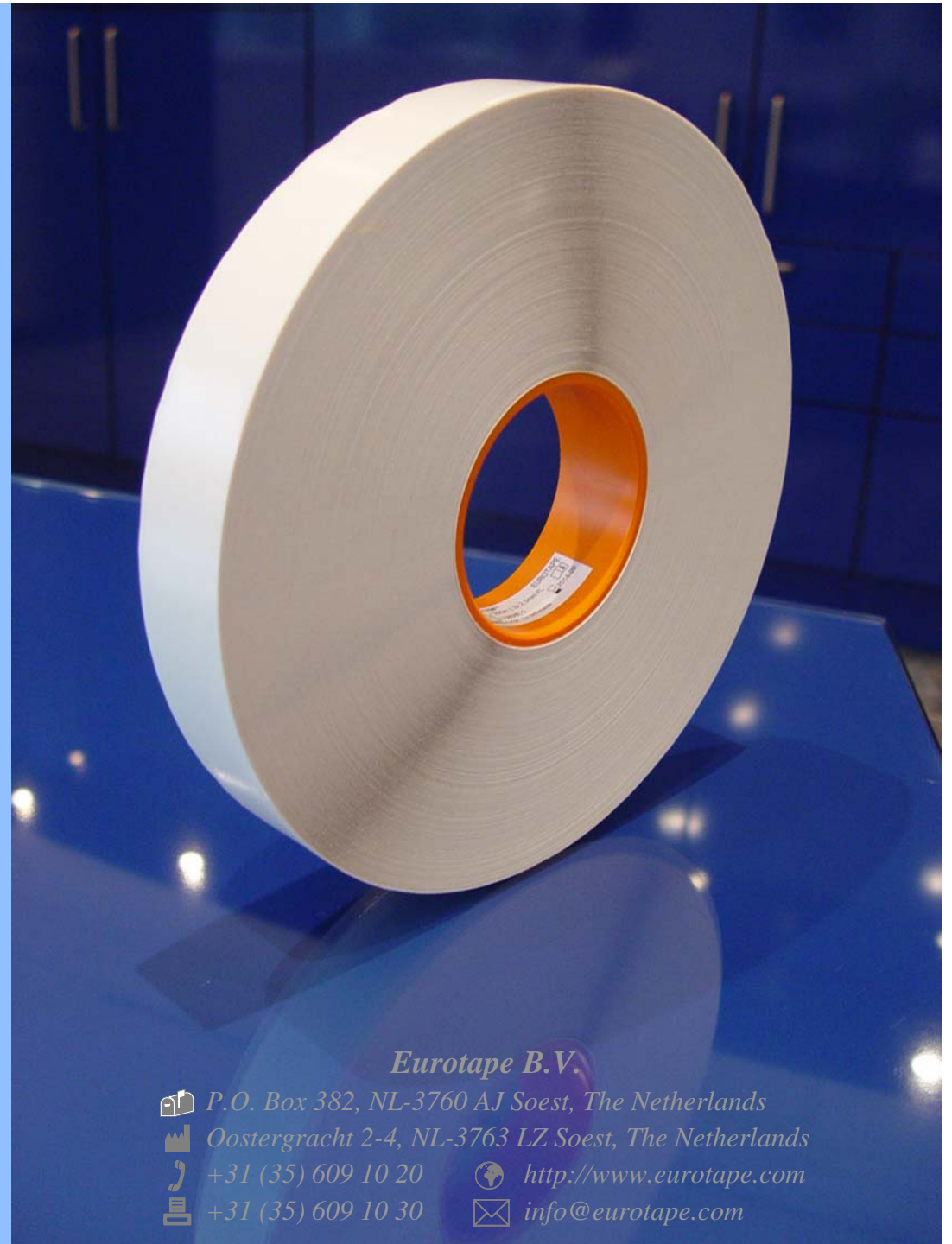

- Certified Management System
- EN ISO 9001
- EN ISO 13485

Eurotape B.V.



P.O. Box 382, NL-3760 AJ Soest, The Netherlands



Oostergracht 2-4, NL-3763 LZ Soest, The Netherlands



+31 (35) 609 10 20



<http://www.eurotape.com>



+31 (35) 609 10 30



info@eurotape.com

Introduction

Intended Use

The OR-tape Med/PV175N™ is intended to be used for applying re-usable drapes to the patient's skin in order to create a sterile and hygienic operating environment.

Product Description


Med/PV175N™ is a double sided tape consisting of an innovative water soluble carrier, both sides of which are coated with a uniform polymer 100% wash soluble adhesive layer. The adhesives are pressure and temperature sensitive and completely wash soluble.

Quality Assurance

All our products are manufactured in accordance with the ISO 9001 and ISO 13485 standards. This ensures that we deliver products and services that consistently meet and exceed customer requirements, allowing for full customer satisfaction.




Our company and products are FDA listed.

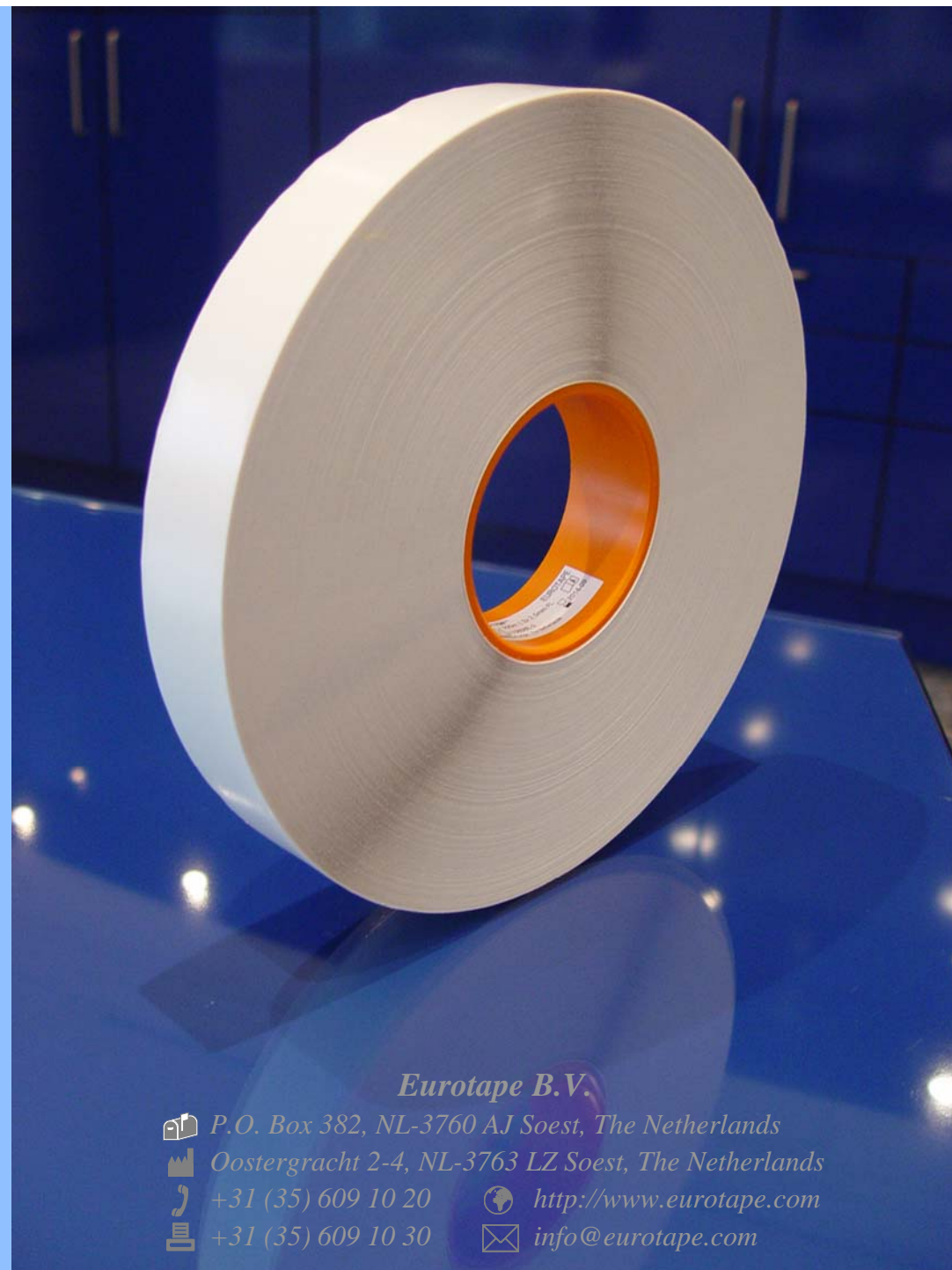We are proud to be the only company in our area of specialization that is entitled to display the CE marking on its products as they comply with the requirements of the Medical Device Directive 2007/47/EEC.




- Certified Management System
- EN ISO 9001
- EN ISO 13485

Eurotape B.V.

 P.O. Box 382, NL-3760 AJ Soest, The Netherlands
 Oostergracht 2-4, NL-3763 LZ Soest, The Netherlands
 +31 (35) 609 10 20  <http://www.eurotape.com>
 +31 (35) 609 10 30  info@eurotape.com

USP's of Eurotape's OR-tape Med/PV175N™

- 100% wash soluble
- 100% biocompatible
- The innovative PVAL carrier does not leave any residues behind after the washing and drying process
- During surgery, it remains unaffected by water and body fluids
- The tape reaches its ultimate adhesive strength on the textile side after steam sterilization
- Custom-made dimensions available upon request

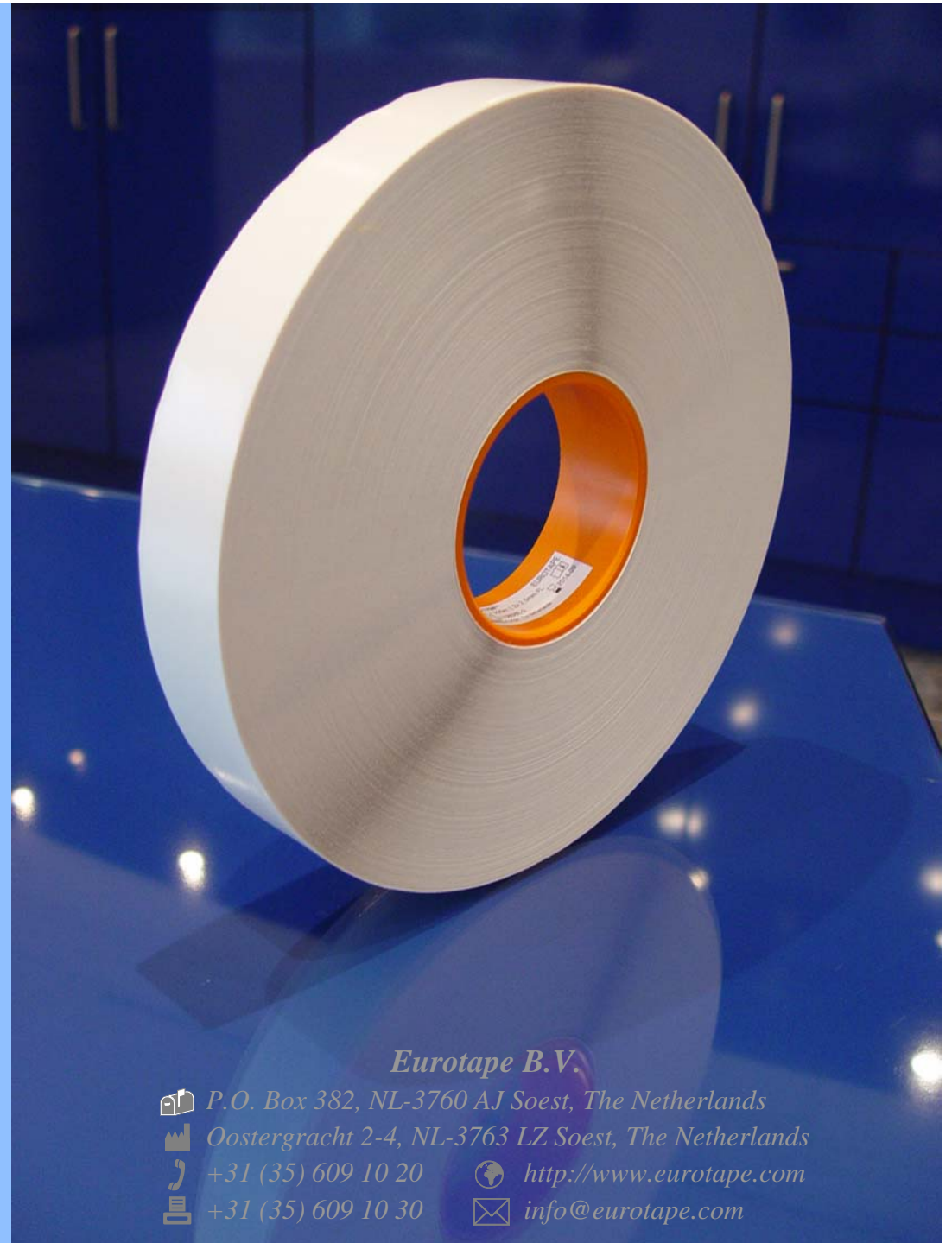

- Certified Management System
- EN ISO 9001
- EN ISO 13485

Eurotape B.V.



P.O. Box 382, NL-3760 AJ Soest, The Netherlands



Oostergracht 2-4, NL-3763 LZ Soest, The Netherlands



+31 (35) 609 10 20



http://www.eurotape.com



+31 (35) 609 10 30



info@eurotape.com